



Contour next link

Wireless Blood Glucose Monitoring System




NO CODING


Uses only CONTOUR[®]NEXT Blood Glucose Test Strips from Bayer



QUICK REFERENCE GUIDE

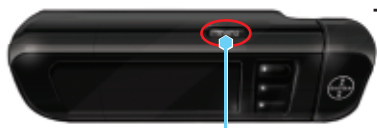

TOP VIEW

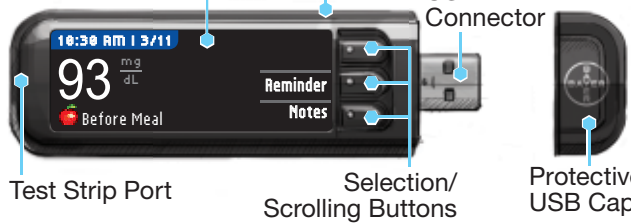

FRONT VIEW

Display Screen

Top Menu Button

USB Connector

Test Strip Port

Selection/Scrolling Buttons

Protective USB Cap

Note: Check the units of measure of your meter. See your meter user guide for complete instructions.

81548021
Rev. 04/13

Initial Setup of the Meter and Connection to Compatible Medtronic MiniMed[®] Device

Charging the Meter

- Remove the cap with the  on it to reveal the USB connector.

The CONTOUR[®]NEXT LINK meter from Bayer has a permanent rechargeable battery. **Before you do anything else**, charge your meter with the wall charger included with your kit or by plugging the USB connector into a computer. Be sure the computer is turned on. Your meter will briefly display “Do Not Test–Charging” and the test strip port light will flash.

You cannot do a blood glucose test while the battery is charging. When charging is finished, the test strip port light will turn off. Unplug your CONTOUR NEXT LINK meter to test.



Setting Up the Meter

- To begin setup, hold down Menu button for several seconds until the meter turns on.
- Choose your language and press **OK**.
- Press **OK** on the next screen to begin the connection process. **NOTE:** Pump refers to your compatible Medtronic MiniMed device.
- STOP. Do NOT press ANY buttons.** This screen displays the meter ID number (in the bottom left corner). This ID is needed to complete the connection process. Put down your meter. Pick up the compatible Medtronic MiniMed device and continue with Step 6.

Connecting Meter to Compatible Medtronic MiniMed Devices

- Press **ACT** to go to the **MAIN MENU** of your compatible Medtronic MiniMed device. Scroll to **Utilities** by pressing the down arrow. Press **ACT**.
- If you are using the Paradigm[®] REAL-Time Insulin Pump, scroll to **Meter Options** from the **UTILITIES MENU** and press **ACT**. If you are using the Guardian[®] REAL-Time Monitor, refer to the user guide for connecting your meter.
- Select **On** and press **ACT**. If you are replacing a meter, go to **Step 9**. If you are adding a new meter, go to **Step 11**.
- If replacing a meter, select **Delete ID** and press **ACT**.
- If you have multiple meter IDs, repeat Step 9 for each one.


- To add a new meter, select **Add ID** and press **ACT**. Refer to Step 5 (or see meter back) to find meter ID. Use the up and down arrow buttons to enter the six character meter ID shown on the meter screen. Press **ACT** after each character entry.



- Select **Review ID** on your compatible Medtronic MiniMed device to confirm that your CONTOUR[®]NEXT LINK meter ID has been correctly entered. If incorrect, go back to Step 9. **NOTE:** For detailed instructions, see your compatible Medtronic MiniMed device user guide.




- Return to your CONTOUR NEXT LINK meter. After you have entered the meter ID into your compatible Medtronic MiniMed device, press **OK**.

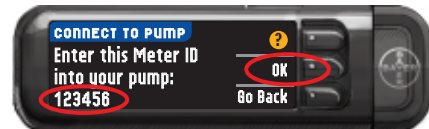

Choosing Meter Send Options

- SEND OPTIONS** controls how your meter will send blood glucose readings to your compatible Medtronic MiniMed device. Choose the **SEND OPTION** that is best for you. **Ask Me** = Meter will ask you before sending results. **Never** = Meter will never send results.

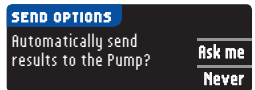

Set Meter Date/Time to Match Medtronic Device

- The date and time on your meter must match the date and time on your Medtronic device. Pick up Medtronic device. Press **ACT** to go to the **MAIN MENU**. Scroll to **Utilities** and press **ACT**. Scroll to **Time/Date** and press **ACT**. Select a time format and press **ACT**. Note the time and date. Return to your meter. Follow instructions on your meter to change the time and date to match your Medtronic device.



Additional Meter Settings

- Choose your **AUTOLOG** settings. AUTOLOG allows you to mark a test as **Before Meal**, **After Meal**, or **No Mark**.
- The meter comes preset with the default blood glucose target ranges as shown. **Accept** the default blood glucose target ranges or **Change** them after talking to your health care professional.
- Follow the display prompts to set your meter **High and Low Limits**. Choose **Accept** to use the default setting (250 mg/dL for the meter high limit and 70 mg/dL for the meter low limit). Choose **Change** if you wish to use different limits.




NOTE: High and Low Limits can only be set outside your selected Target range(s). Targets can only be set inside your selected High and Low Limit(s). For additional information, you may contact Customer Service at 800-348-8100, available 24 hours a day, 7 days a week.

SETUP IS NOW COMPLETE!

Materials for Testing

 CONTOUR®NEXT
Test Strip from Bayer

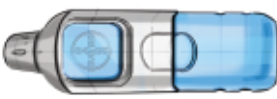 MICROLET®2
Lancing Device
from Bayer

 MICROLET®
Lancet from Bayer

 CONTOUR®NEXT LINK
Blood Glucose Meter
from Bayer

 WARNING:


- The meter, lancing device and lancets are for single patient use. Do not share them with anyone including other family members. Do not use on multiple persons.^{1,2}
- All parts of the kit are considered biohazardous and can potentially transmit infectious diseases, even after you have performed cleaning and disinfection.
- Always wash your hands well with soap and water before and after testing, handling the meter, lancing device, or test strips.


For complete instructions on cleaning and disinfecting your meter and lancing device, please see your meter user guide.

Testing Your Blood

Obtaining the Blood Drop

- 1 Insert test strip and wait for **APPLY BLOOD** message.

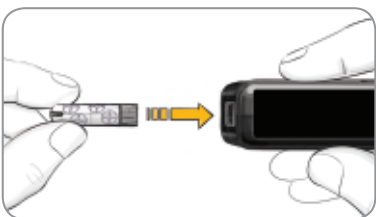
- 2 Remove lancing device endcap. Insert lancet into lancing device and twist off protective cap.

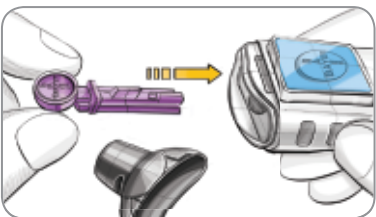
- 3 Re-cap and dial to desired lancing depth.


- 4 Obtain your blood drop.



Testing the Blood Drop

- 5 Touch the tip of the test strip to the drop of blood. If the first blood drop is not enough, the meter will beep twice and display **STRIP UNDERFILLED**. You have about 30 seconds to apply more blood to the same strip.


- 6 If AUTOLOG is on, you can choose to record your result as:  **Before Meal**,  **After Meal** or **No Mark** for the test result.



Sending Results to Your Compatible Medtronic MiniMed Device

- 7 If **Ask Me** was selected in **Send Options** during **Setup**, results can be sent to your compatible Medtronic MiniMed device. Press **OK** to move to the next screen.




Reminder and Notes

- 8 Select **Reminder** to set a reminder for your next test. Select **Notes** to add more information to the result. Remove test strip.



Turn Meter Off

- 9 Remove the test strip and the meter will turn off in 3 minutes. Or hold down the Menu button to turn off the meter.

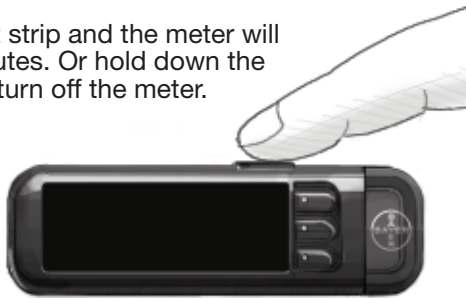

See your CONTOUR®NEXT LINK user guide, MICROLET®2 insert and compatible Medtronic MiniMed device user guide for full instructions.

Call **800-348-8100** (available 24 hours a day, 7 days a week) with any questions

or visit www.bayerdiabetes.com/us

For In Vitro Diagnostic Use

References

- ¹ FDA Public Health Notification: Use of Fingerstick Devices on More than One Person Poses Risk for Transmitting Bloodborne Pathogens: Initial Communication (2010). <http://www.fda.gov/MedicalDevices/Safety/AlertsandNotices/ucm224025.htm>
- ² CDC Clinical Reminder: Use of Fingerstick Devices on More than One Person Poses Risk for Transmitting Bloodborne Pathogens (2010). <http://www.cdc.gov/injectionsafety/Fingerstick-DevicesBGM.html>

 **Bayer HealthCare**

Bayer HealthCare LLC
Mishawaka, IN 46544 USA

Bayer, the Bayer Cross, CONTOUR, MICROLET and the No Coding logo are registered trademarks of Bayer. All other trademarks are the property of their respective owners.

© 2013 Bayer. All rights reserved.

 **Medtronic**